

## Central Lancashire Online Knowledge (CLoK)

Title	Collection of Touch DNA from Rotten Banana Skin
Type	Article
URL	<a href="https://clock.uclan.ac.uk/37762/">https://clock.uclan.ac.uk/37762/</a>
DOI	10.23880/ijfsc-16000204
Date	2020
Citation	Alketbi, Salem Khalifa (2020) Collection of Touch DNA from Rotten Banana Skin. International Journal of Forensic Sciences, 5 (4). ISSN 2573-1734
Creators	Alketbi, Salem Khalifa

It is advisable to refer to the publisher's version if you intend to cite from the work. 10.23880/ijfsc-16000204

For information about Research at UCLan please go to <http://www.uclan.ac.uk/research/>

All outputs in CLoK are protected by Intellectual Property Rights law, including Copyright law. Copyright, IPR and Moral Rights for the works on this site are retained by the individual authors and/or other copyright owners. Terms and conditions for use of this material are defined in the <http://clock.uclan.ac.uk/policies/>



# Collection of Touch DNA from Rotten Banana Skin

Salem K Alketbi<sup>1,2\*</sup>

<sup>1</sup>University of Central Lancashire, UK

<sup>2</sup>General Department of Forensic Science and Criminology, UAE

**\*Corresponding author:** Salem K Alketbi, University of Central Lancashire, Preston, UK, General Department of Forensic Science and Criminology, Dubai Police, UAE, Tel: 00447774141205; Email: alkitbe.11@hotmail.com

## Research Article

Volume 5 Issue 4

Received Date: October 08, 2020

Published Date: November 02, 2020

DOI: 10.23880/ijfsc-16000204

## Abstract

Touch or trace DNA analysis can be used to connect suspects to crime scenes, especially when other types of biological evidence are unavailable. However, Touch DNA profiling is a challenging process and many variables affect the success of obtaining a good quality DNA profile, such as surface type and the time between deposition and recovery. This study aimed to investigate the influence of time between deposition and recovery of Touch DNA from banana skin surface at room temperature. The results showed that the amount of recovered DNA from the banana skin was not affected by time (over a week) at room temperature and a full DNA profile was generated even when the banana skin was rotten.

**Keywords:** Forensic DNA; Touch DNA; DNA recovery; Copan 150C Cotton swab; PrepFiler Express BTA™; QIAamp® DNA Investigator; Quantifiler™ Human DNA Quantification Kit; QuantStudio; GlobalFiler™; PCR Amplification Kit

## Introduction

Touch or trace DNA analysis is a significant tool that can be used to connect suspects to crime scenes, especially when other types of biological evidence cannot be found. The process of Touch DNA profiling has evolved since it was first reported in 1997 [1] with the improvement of DNA profiling technology. Indeed, it has become an important routine of the forensic laboratory workload and unlocked new possibilities that led to the collection of DNA from a wide range of surfaces and items, such as tools, knives, clothing, firearms, etc. [2,3].

Touch DNA profiling is challenging, and many variables can affect the success of obtaining a good quality DNA profile, such as surface type and the time between deposition and recovery [4-6]. Previous studies have mainly focused on improving DNA recovery from body fluids [7,8], thus there is a gap in the research regarding Touch DNA [3].

During the examination of numerous home burglary crime scenes, while working for the General Department

of Forensic Science and Criminology in Dubai Police, it was noted that some criminals consume food from the home, such as banana, and then leave the banana skin at the crime scene. Therefore, this study aimed to investigate the influence of time between deposition and recovery of Touch DNA from banana skin surface at room temperature.


## Materials and Methods

### Experimental Setup and Deposition

To assess the effect of time on the collection of Touch DNA from banana skin, it was stored at room temperature (RT) in moderate humidity (20 to 25°C/50%) (n=15) over five time periods (Zero hours, 3 hours, 12 hours, 24 hours and one week). This was done to test the possibility of collecting Touch DNA from rotten banana skin (Figure 1).

Before storing the items, surfaces were cleaned with viricidal disinfectant (2% virkon) and ultraviolet radiation (UV) for 15 min. The participant was asked to wash their

hands with antibacterial soap and refrain from undertaking any activity for 10 minutes, then charge the fingers of both hands with eccrine sweat by touching behind their ears or forehead to load them with enough DNA. The participant was then asked to touch the surfaces using their index, middle, and ring fingers of both hands separately to deposit the DNA by applying medium pressure on a 5 x 7 cm surface area for 1 minute. The same procedure was repeated for each deposition.

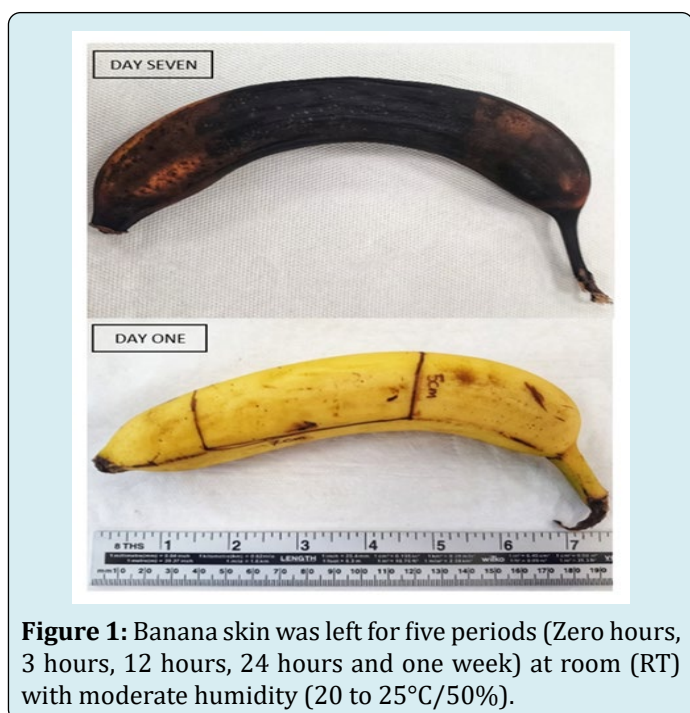

**Figure 1:** Banana skin was left for five periods (Zero hours, 3 hours, 12 hours, 24 hours and one week) at room (RT) with moderate humidity (20 to 25°C/50%).

### DNA Recovery and Extraction

Samples were collected using a Copan 150C Cotton


swab (Copan, Brescia, Italy) moistened with 100  $\mu\text{L}$  of sterile distilled water applied using a plastic spray bottle technique (developed in Dubai police forensic DNA lab; every single spray contains approximately 50  $\mu\text{L}$ ). Full swab heads were used for extraction immediately after collection using the PrepFiler Express BTA™ kit with AutoMate Express Forensic DNA Extraction System according to the manufacturer's instructions (Thermo Fisher Scientific) and the final extracted sample elution was 50  $\mu\text{L}$ .

### DNA Quantification, Amplification and Analysis


Extracted samples were quantified using the Quantifiler® Human DNA Quantification Kit, QuantStudio 5 Real-Time PCR (qPCR) and HID Real-Time PCR analysis software v1.3 according to the manufacturer's instructions (Thermo Fisher Scientific). Amplification was performed using the GlobalFiler™ PCR Amplification Kit (Thermo Fisher Scientific) and a 30-cycle protocol. The data were analysed using GeneMapper® ID-X Software Version 1.2 (Thermo Fisher Scientific) and statistical analysis was performed on the tested variables using RStudio and factorial analysis of variance (ANOVA). In ANOVA, the p-value is derived from the F-distribution which is different for every pair of degrees of freedom (df) values ( $F \text{ value} = \text{variance of the variables means (Mean Square Between)} / \text{mean of the within variables variances (Mean Squared Error)}$ ). Blanks were taken from the surfaces after sterilisation, and negative controls for the collection and extraction methods, all of which proved negative for DNA when quantified.

### Results and Discussion

The amount of recovered DNA from the banana skin was not affected by time (over a week) at RT ( $F_{4,5} = 0.33$ ,  $p > 0.05$ , mean = 0.03 ng/ $\mu\text{L}$ ) (Figure 2), and a full DNA profile was generated even when the banana skin was rotten (Figure 3).

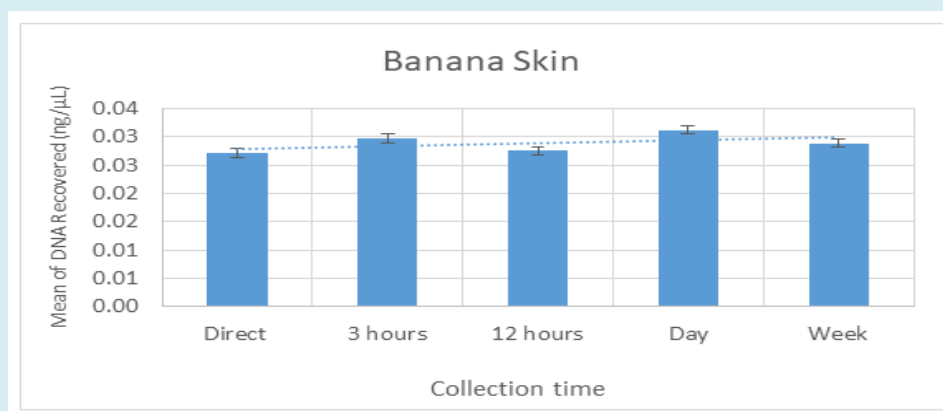

**Figure 2:** The mean of DNA collected from banana skin, at five periods (Zero hours (direct collection after deposition), 3 hours, 12 hours, 24 hours and one week) ( $n=15$  - three replicates). Error bars represent 95% confidence intervals.

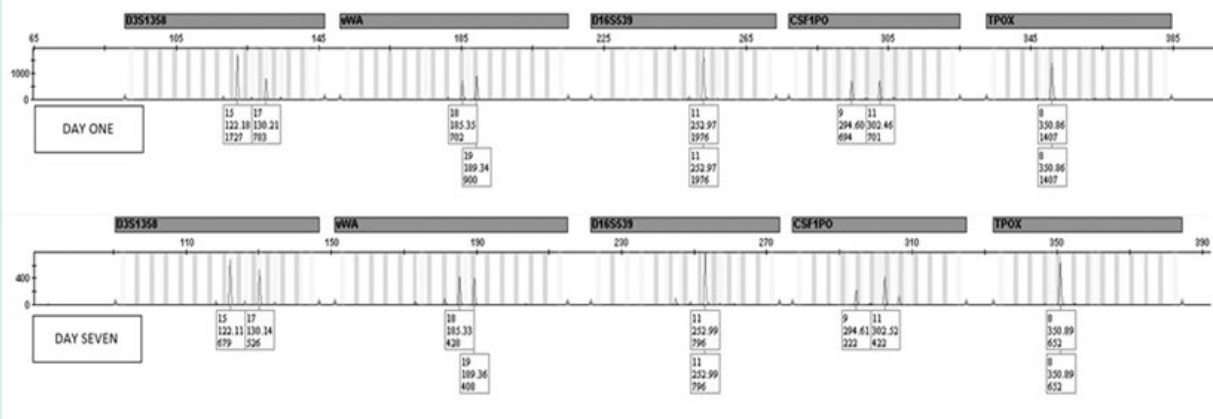

**Figure 3:** Electropherograms of samples collected from normal banana skin (day one) and rotten banana skin (day seven). The profiles show the difference in peaks height at 5 autosomal STR loci (D3S1358, vWA, D16S539, CSF1PO, TPOX).

## Conclusion

In conclusion, this experiment showed that collection of Touch DNA from fresh or rotten banana skin is possible if the banana skin surface was found indoors at room temperature with moderate humidity (20 to 25°C/50%) within one week of deposition. It is important to consider other elements that could influence the Touch DNA profiling process from surfaces, such as environmental conditions and contamination.

**Conflict of Interest:** None.

## Acknowledgement

This study was approved by the General Department of Forensic Science and Criminology in Dubai Police and ethical approval was granted by the School of Forensic and Applied Sciences, and the University of Central Lancashire's Research Ethics Committee (ref. no. STEMH 912). Many thanks to COPAN DIAGNOSTICS INC. for supporting this experiment with free swabs, and to ThermoFisher Scientific™ for the discounts on their products.

## References

1. Van Oorschot RAH, Jones MK (1997) DNA fingerprints from fingerprints. *Nature* 387: 767.
2. Polley D, Mickiewicz P, Vaughn M, Miller T, Warburton R, et al. (2006) An investigation of DNA recovery from firearms and cartridge cases. *J Canadian Soc Forensic Sci* 39(4): 217-228.
3. Alketbi SK (2018) The affecting factors of Touch DNA. *Journal of Forensic Research* 9(3): 424.
4. Bille TW, Cromartie C, Farr M (2009) Effects of cyanoacrylate fuming, time after recovery, and location of biological material on the recovery and analysis of DNA from post, blast pipe bomb fragments. *J Forensic Sci* 54(5): 1059-1067.
5. Alketbi SK, Goodwin W (2019) The effect of surface type, collection and extraction methods on touch DNA. *Forensic Science International. Genetics Supplement Series* 7(1): 704-706.
6. Alketbi SK, Goodwin W (2019) The effect of time and environmental conditions on Touch DNA. *Forensic Science International. Genetics Supplement Series* 7(1): 701-703.
7. Verdon T, Mitchell R, Oorschot R (2014) Swabs as DNA Collection devices for sampling different biological materials from different substrates. *Journal of Forensic Sciences* 59(4): 1080-1089.
8. Plaza D, Mealy J, Lane J, Parsons M, Bathrick A, et al. (2016) Nondestructive biological evidence collection with alternative swabs and adhesive lifters. *Journal of Forensic Sciences* 61(2): 485-488.

